Compression Garments by NordiCare

# Kim

art no 51107

POST SURGICAL FACIAL BANDAGE

2021-11-11

Available colors: 

 This packaging can be recycled as paper and the bag as plastic.

 NordiCare Ortopedi & Rehab AB  
Sörsvägen 1,  
SE-263 62 Viken  
Phone +46 42 35 22 20,  
Fax +46 42 35 22 21  
info@nordicare.se,  
www.nordicare.se

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**EN** Facial Bandage in comfortable material with seamless edges. Adjustable fit and compression. The material is possible to trim/cut to allow optimum comfort and fit.

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**EN** Possible to trim and cut the material and straps for optimal fit. For example, trim around the cheeks, chin and neck, or cut holes for the ears.

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**EN** At facelift and/or liposuction of the chin/neck, ear surgery or when the steady compression under the chin, around the neck, cheeks, hairline or ears is desired. **Contraindications:** None known.

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
**EN** Machine wash at 40°C (104°F), can if necessary be washed at higher temperatures, up to 60°C (140°F). The product will last longer if washed at a lower temperature. Dry flat. Do not tumble dry. Close the Velcro fastening before washing, use a washing net. In case of blood stains, soak or wash in cold water. Do not use bleach or fabric softener as this can impair the function of the material.

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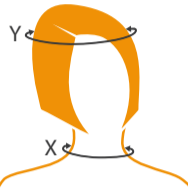**Kaula:** Kaulan levein kohta.

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**EN** Use the size guide to decide the correct size. The sizes are flexible and can be adjusted and trimmed for desired compression and optimal fit.

**Neck:** The widest point in the middle of the neck.

**Head:** Widest point over the forehead.

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## Materials

Polyamide 76%, Elastane 14%, Polyurethane 10%. Latex free.

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**EN** If a stinging sensation, rash or other changes in skin color occur, discontinue use and contact your doctor or other medically trained person. Only medically trained persons are to prescribe and recommend the use of compression garments. For best results, follow instructions for use. Serious incidents shall be reported to NordiCare and the competent authority.

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Please consult our product catalogue for more information. For further questions, contact your doctor, health practitioner, medical equivalent or NordiCare: +46 42-35 22 20 / DK 7070-1907, info@nordicare.se.